

# Sampling Bottle 3



- For Quantitative Fecal Immunochemical Test (FIT)
- For Quantitative Fecal Calprotectin Test (FCa)

## Ease, Accuracy and Hygiene all in one compact body

### Flat Surface

Allows for easy writing of patient information

### Sample Buffer

Stabilizes the sample for an extended period

### Dual Aluminum Seal

Prevents accidental leakage for improved hygiene

### Integrated Filter

Removes fecal matter from sample solution through the squeeze filtration process

### Scraper

Ensures accurate quantitative results by removing excess feces in the chamber

### Probe Groove

Provides smooth and precise sampling of 10mg feces

### Hollow Cap

Easy to open and close, preventing sample leakage and contamination hazards



## Unique Squeeze Filtration Process

Only filtered sample solution is dispensed for analysis.

This mechanism minimizes carry-over and analyzer maintenance, as well as helps achieve high reproducibility.

# For Fecal Collection and Sample Return



## Collection Paper

- Printed sampling instruction
- Flushable paper
- Available in different designs




## Special Envelope


- Integrated absorbent pad
- Available for both 1 and 2 day sampling methods

# Stability in the Buffer



## Hemoglobin

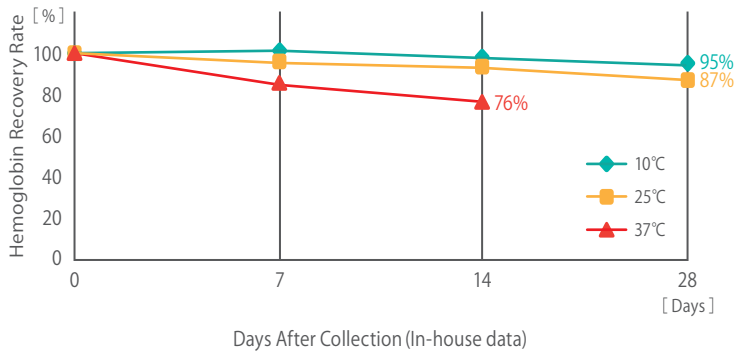

## Calprotectin

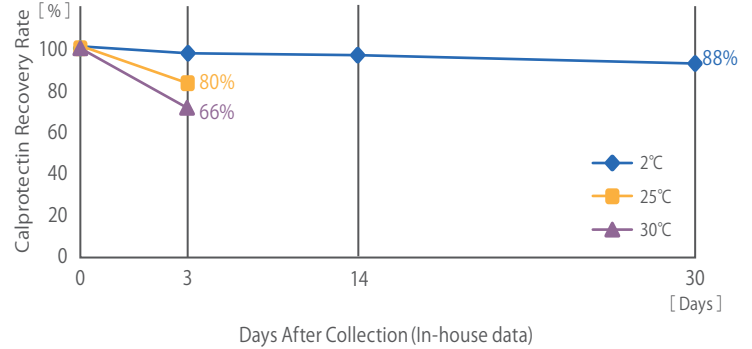

# How to Collect Sample

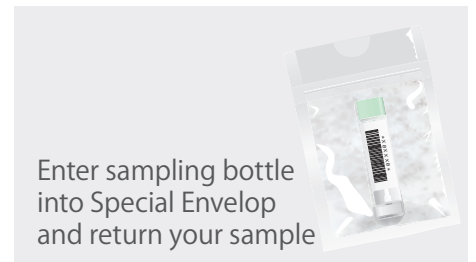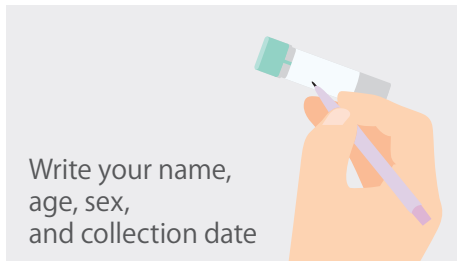

Ref No.	Product name	Contents	Storage
V-PZ25	OC-Auto Sampling Bottle 3	100 bottles	1-30°C
V-PZ26	OC-Auto Sampling Bottle 3 without barcode	100 bottles	1-30°C